

pulmoNeb

Small & robust compressor nebuliser



The 3655ix PulmoNeb compressor nebuliser from DeVilbiss is small, robust and simple to use.

Its compact design includes a wrap around power cord for ease of storage, built in nebuliser holder, replaceable filter and carry handle. The air outlet port located on the top of the unit makes it easy for the user to attach and remove the tubing. The PulmoNeb is suitable for both adults and children and can be used with a disposable or reusable nebuliser such as the optional DeVilbiss Neb 800.

Features & Benefits

- Small and robust design
- Easy to operate and clean
- Available with both reusable and disposable nebulisers
- Easily accessible single piece replaceable filter

- Metal mechanical parts
- Wrap around power cord storage feature

Neb 800 – optional reusable nebuliser

- 75 % of the medication is deposited in the lungs with an MMAD of 3.6 – 4 microns
- High aerosol output to reduce treatment time
- Spill proof design
- Long term usage up to one year
- Nebuliser can be cleaned in the dishwasher
- Mouthpiece can be replaced with a standard mask
- Less medication waste due to effective treatment and residual volume < 0.5 ml
- Gives effective treatment and results, even at a 45° angle



SPECIFICATIONS

Dimensions (H x W x D)

19.1 x 10.2 x 17.3 cm

Weight 1.86 kg

Electrical requirements 230 V / 50 Hz

Maximum pressure
35 psig or greater

Compressor free air flow
9 l/min or greater

Neb 800:

Pressure Rate Nebuliser 1 - 1.4 bar

Nebuliser Litre Flow 5 - 8 l/min.

Nebuliser Capacity 6 ccm

Nebuliser Rate 0.3 ml/min.

Particle Size / MMAD


75% < 5 µm / 3.6 - 4.0 µm



PulmoNeb compressor nebuliser



PulmoNeb with reusable Neb800

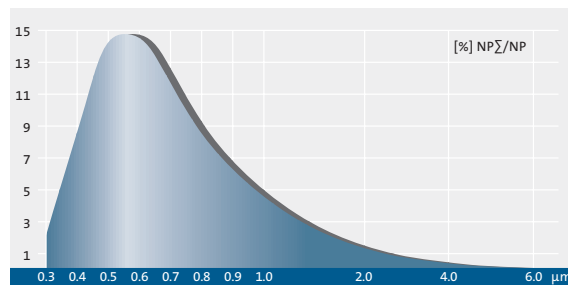

Particle Size



Neb 800 – optional reusable nebuliser