

50A


STERYBOX®

Professional UVC Airborne Disinfection Device

Made in Italy

Ultraviolet germicidal irradiation (UVGI) has a long history of studies showing it is safe and highly-effective in air sterilization. Lots of clinical investigations conclude UVGI is good at preventing transmission of a variety of infectious diseases, evidenced by the very significant drop of infection rate when UVGI was in use.

Recently, the Italian-made Sterybox®, which manufacturer is **an expert in UVC air sterilizer for over 25 years**, was introduced to Hong Kong. It is the only device in the market with official microbiological decontamination results and dedicated scientific documentation performed by independent laboratories recognized by Government Health Institutions. Installation of thousands of units per year in hospitals and operation theatres is a guarantee itself.

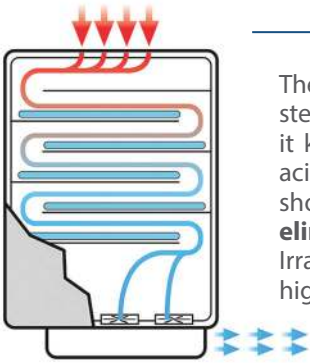

The Sterybox® consists of FIVE sets of UV lamps in a sealed stainless steel box. Using the shortwave UV-C (at 254 nm) emitted by the lamps, it kills and inactivates microbes by damaging their deoxyribonucleic acid (DNA). EU accredited microbiology lab results on Sterybox® showed **99.99% of airborne bacteria, viruses and micro spores were eliminated**. Its patented UV-CHRIS design (UltraViolet High Reflective Irradiation Coil System) delivers a UV concentration up to four times higher than the standard linear UV lamps, allowing air passing through a long 2.8-metre UV irradiation path before releasing out back to the room fresh and cleaned.

Sterybox® is safely sealed by means of various panels. The germicidal irradiation remains only inside the unit and the operators can work absolutely safe in the room when Sterybox® is on. There is also a safety anti-opening switch control that switches off the UV light in the event that the box is accidentally opened. Therefore, it eliminates any potential danger from long UV-C exposure.

No risky handling or replacement of contaminated filters and quiet operation with powerful internal silencer, Sterybox® is a perfect convenient choice for many professional practices.



Ambulance

Clinic

Dental Office

Sterybox® is exclusively available at The HomeCare Medical Ltd now.

For detailed information, please call 2402-2188 or visit our website <https://homecare-medical.com/en/other/respiratory/sterybox/>

